

Integra™ Jarit®

Sterilization Containers

Limit uncertainty by protecting
your instrument investment



INTEGRA™
LIMIT UNCERTAINTY



Advantages of Integra™ Jarit® Containers over Wrap

- Clinically advantageous method to effectively sterilize, store and transport surgical instruments
- Represents a cost-saving opportunity
- Reduces labor and provides improved service to both SPD and the OR
- Environmentally – friendly
- ID tags offer identification of contents and destination
- Easy to lift and carry
- Improves set preparation times

Summary of Integra Jarit Sterilization Container Product Features

Silicone gaskets within retention plates and around container perimeter ensure a secure seal of the container, lid and filter.

Snap on filter plate with a simple audible click and lock mechanism for quick, secure and easy change of disposable filters.

Round disposable filters with centering hole to assist in proper alignment.

Lid and Bottom manufactured from anodized aluminum alloy. Easy to maintain, corrosion resistant, lightweight yet durable, conducts heat for material drying and prevents passivation between stainless steel instruments and the container itself.

4 ID Tag Windows (2 front & 2 back)

ID tags available in 6 colors and can be engraved with up to 18 characters.


Low profile latch easily and securely closes lid over container bottom.

Tamper Evident Lock – identifies sterility breach.

Fold-down handles extend and stop 90° for secure grasp, aids lifting and avoids pinching.

Corner design allows for same size container stacking in storage.

2 Indicator Card Windows (1 front & 1 back). Indicator card has visible bar for clear indication that steam or gas sterilization has been achieved.

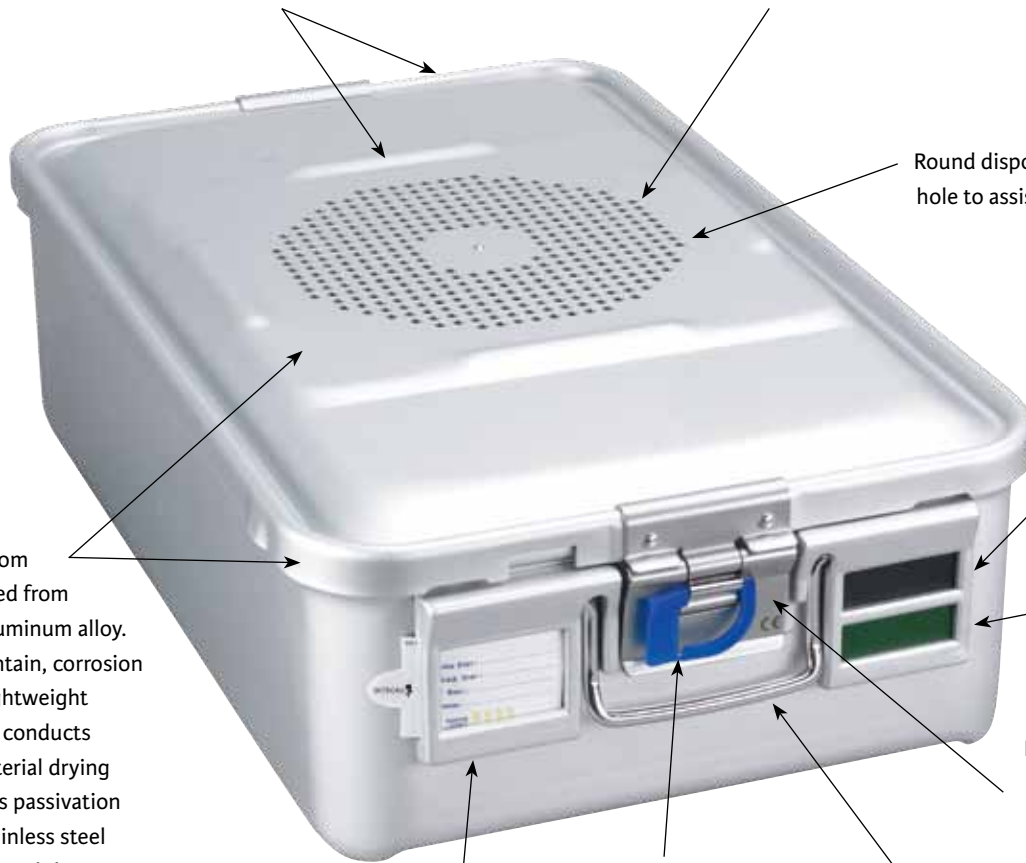

*Baskets are constructed with solid wire mesh to allow for maximum ultrasonic cleaning; handles can be positioned to assist with secure and comfortable grasp for removal under sterile conditions.

One-Half Size (11-1/4" x 11") Solid Container Bottom with Lid

- Vacuum Sterilizer Containers - Solid/Perforated Bottom
- Gravity Sterilizer Containers - Perforated Bottom
- Ethylene Oxide Sterilizer Containers - Solid/Perforated Bottom

Catalog Number	Height	Description
730-100	4"	One-half size, solid container bottom with lid, silver
730-101	5 1/3"	One-half size, solid container bottom with lid, silver
730-102	6"	One-half size, solid container bottom with lid, silver
730-103	8"	One-half size, solid container bottom with lid, silver

Catalog Number	Height	Description
730-110	4"	One-half size, perforated container bottom with lid, silver
730-111	5 1/3"	One-half size, perforated container bottom with lid, silver
730-112	6"	One-half size, perforated container bottom with lid, silver
730-113	8"	One-half size, perforated container bottom with lid, silver



Note: Complete container assembly only available in silver.
See below for other color options with lids and bottoms sold separately.

One-Half Size (11-1/4" x 11") Colored Lids (bottoms sold separately)

Includes Filter Retention Plate and Latches

Catalog Number	Description
730-120	Jarit colored, lid container, one-half size, silver
730-121	Jarit colored, lid container one-half size, blue
730-122	Jarit colored, lid container one-half size, red
730-123	Jarit colored, lid container, one-half size, gold
730-124	Jarit colored, lid container-half size, green
730-125	Jarit colored, lid container, one-half size, black



One-Half Size Container Bottoms (lids sold separately)

Catalog Number	Height	Description
730-100B	4"	One-half size, solid container bottom only
730-101B	5 1/3"	One-half size, solid container bottom only
730-102B	6"	One-half size, solid container bottom only
730-103B	8"	One-half size, solid container bottom only



Catalog Number	Height	Description
730-110B	4"	One-half size, perforated container bottom only
730-111B	5 1/3"	One-half size, perforated container bottom only
730-112B	6"	One-half size, perforated container bottom only
730-113B	8"	One-half size, perforated container bottom only



One-Half Size Wire/Perforated Baskets

Wire Baskets

Catalog Number	Height	Description
740-100	2"	One-half size wire baskets
740-101	3"	One-half size wire baskets
740-102	4"	One-half size wire baskets

Feet help to raise basket over retention plate and are ideal for use with gravity perforated bottom containers. Feet are removable.
Note: Plastic feet for baskets, PKG/1 (Spare Part: 780-FT-100)



Perforated Baskets

Catalog Number	Height	Description
750-100	2"	One-half size perforated baskets
750-101	3"	One-half size perforated baskets
750-102	4"	One-half size perforated baskets

Feet help to raise basket over retention plate and are ideal for use with gravity perforated bottom containers. Feet are removable.
Note: Plastic feet for baskets, PKG/1 (Spare Part: 780-FT-100)



Three-Quarter Size (18-1/4" x 11") Solid Container Bottom with Lid

- Vacuum Sterilizer Containers - Solid/Perforated Bottom
- Gravity Sterilizer Containers - Perforated Bottom
- Ethylene Oxide Sterilizer Containers - Solid/Perforated Bottom

Catalog Number	Height	Description
730-200	4"	Three-quarter size, solid container bottom with lid, silver
730-201	5-1/3"	Three-quarter size, solid container bottom with lid, silver
730-202	6"	Three-quarter size, solid container bottom with lid, silver

Catalog Number	Height	Description
730-210	4"	Three-quarter size, perforated container bottom with lid, silver
730-211	5-1/3"	Three-quarter size, perforated container bottom with lid, silver
730-212	6"	Three-quarter size, perforated container bottom with lid, silver



Note: Complete container assembly only available in silver.
See below for other color options with lids and bottoms sold separately.

Three-Quarter Size (18-1/4" x 11") Colored Lid (bottoms sold separately)

Includes Filter Retention Plate and Latches

Catalog Number	Description
730-220	Jarit colored, lid container, three-quarter size, silver
730-221	Jarit colored, lid container, three-quarter size, blue
730-222	Jarit colored, lid container, three-quarter size, red
730-223	Jarit colored, lid container, three-quarter size, gold
730-224	Jarit colored, lid container, three-quarter size, green
730-225	Jarit colored, lid container, three-quarter size, black



Three-quarter Container Bottoms (lids sold separately)

Catalog Number	Height	Description
730-200B	4"	Three-quarter size, solid container bottom only
730-201B	5-1/3"	Three-quarter size, solid container bottom only
730-202B	6"	Three-quarter size, solid container bottom only



Catalog Number	Height	Description
730-210B	4"	Three-quarter, perforated container bottom only
730-211B	5 1/3"	Three-quarter, perforated container bottom only
730-212B	6"	Three-quarter, perforated container bottom only



Three-Quarter Size Wire/Perforated Baskets

Wire Baskets

Catalog Number	Height	Description
740-200	2"	Three-quarter size wire baskets
740-201	3"	Three-quarter size wire baskets
740-202	4"	Three-quarter size wire baskets

Feet help to raise basket over retention plate and are ideal for use with gravity perforated bottom containers. Feet are removable.

Note: Plastic feet for baskets, PKG/1 (Spare Part: 780-FT-100)




Perforated Baskets

Catalog Number	Height	Description
750-200	2"	Three-quarter size perforated baskets
750-201	3"	Three-quarter size perforated baskets
750-202	4"	Three-quarter size perforated baskets

Feet help to raise basket over retention plate and are ideal for use with gravity perforated bottom containers. Feet are removable.

Note: Plastic feet for baskets, PKG/1 (Spare Part: 780-FT-100)

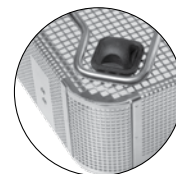

Full-Size (23" x 11") Solid Container Bottom with Lid

- Vacuum Sterilizer Containers - Solid/Perforated Bottom
- Gravity Sterilizer Containers - Perforated Bottom
- Ethylene Oxide Sterilizer Containers - Solid/Perforated Bottom

Catalog Number	Height	Description
730-300	4"	Full-size, solid container bottom with lid, silver
730-301	5-1/3"	Full-size, solid container bottom with lid, silver
730-302	6"	Full-size, solid container bottom with lid, silver
730-303	8"	Full-size, solid container bottom with lid, silver

Catalog Number	Height	Description
730-310	4"	Full-size, perforated container, silver
730-311	5-1/3"	Full-size, perforated container, silver
730-312	6"	Full-size, perforated container,, silver
730-313	8"	Full-size, perforated container,, silver

Note: Complete container assembly only available in silver.
See below for other color options with lids and bottoms sold separately.



Full-Size (23" x 11") Colored Lids (bottoms sold separately)

Includes Filter Retention Plate and Latches

Catalog Number	Description
730-320	Jarit colored, lid container, full-size, silver
730-321	Jarit colored, lid container, full-size, blue
730-322	Jarit colored, lid container, full-size, red
730-323	Jarit colored, lid container, full-size, gold
730-324	Jarit colored, lid container, full-size, green
730-325	Jarit colored, lid container, full-size, black



Container Bottoms (lids sold separately)

Catalog Number	Height	Description
730-300B	4"	Full-size, solid container bottom only
730-301B	5-1/3"	Full-size, solid container bottom only
730-302B	6"	Full-size, solid container bottom only
730-303B	8"	Full-size, solid container bottom only



Catalog Number	Height	Description
730-310B	4"	Full-size, perforated container bottom only
730-311B	5-1/3"	Full-size, perforated container bottom only
730-312B	6"	Full-size, perforated container bottom only
730-313B	8"	Full-size, perforated container bottom only



Full-Size Wire/Perforated Baskets

Wire Baskets

Catalog Number	Height	Description
740-300	2"	Full-size size wire baskets
740-301	3"	Full-size size wire baskets
740-302	4"	Full-size size wire baskets

Feet help to raise basket over retention plate and are ideal for use with gravity perforated bottom containers. Feet are removable.

Note: Plastic feet for baskets, PKG/1 (Spare Part: 780-FT-100)




Perforated Baskets

Catalog Number	Height	Description
750-300	2"	Full-size perforated baskets
750-301	3"	Full-size perforated baskets
750-302	4"	Full-size perforated baskets

Feet help to raise basket over retention plate and are ideal for use with gravity perforated bottom containers. Feet are removable.

Note: Plastic feet for baskets, PKG/1 (Spare Part: 780-FT-100)

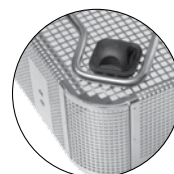

Instrument Stringers Container Accessories

Catalog Number	Description	Size
117-240	Weinstein Instrument Holders	6"
117-242	Weinstein Instrument Holders	8"
117-244	Weinstein Instrument Holders	10"
117-246	Weinstein Instrument Holders	12"
117-247	Weinstein Instrument Holders	14"
117-248	Weinstein Instrument Holders	16"
117-249	Weinstein Instrument Holders	18"
117-250	JARIT Instrument Holders	6"
117-252	JARIT Instrument Holders	8"
117-254	JARIT Instrument Holders	10"
117-256	JARIT Instrument Holders	12"
117-258	JARIT Instrument Holders	14"
117-260	JARIT Instrument Holders	16"
117-262	JARIT Instrument Holders	18"
117-270	Instrument Stringers	8"
117-272	Instrument Stringers	10"
117-274	Instrument Stringers	12"
117-276	Instrument Stringers	14"
117-278	Instrument Stringers	16"
117-280	Instrument Stringers	18"



Note: For more instruments stringers please refer to the Jarit® catalog.

117-240 to 117-249

117-250 to 117-262

117-270 to 117-280

Identification Tags

Catalog Number	Description
775-100	ID Tag, silver
775-101	ID Tag, blue
775-102	ID Tag, red
775-103	ID Tag, gold
775-104	ID Tag, green
775-105	ID Tag, black

Each tag allows for 2 lines of custom text max, 18 characters each.



Small Snap-on Basket Tags (with key ring)

Catalog Number	Description	Size	Size
775-120	Silver	2 3/8" x 1/2"	2 3/8" x 1 1/8"
775-121	Gold	2 3/8" x 1/2"	2 3/8" x 1 1/8"
775-122	Red	2 3/8" x 1/2"	2 3/8" x 1 1/8"
775-123	Blue	2 3/8" x 1/2"	2 3/8" x 1 1/8"
775-124	Green	2 3/8" x 1/2"	2 3/8" x 1 1/8"
775-125	Black	2 3/8" x 1/2"	2 3/8" x 1 1/8"


Each tag allows for 3 lines of custom text max, 32 characters each.



Large Snap-on Basket Tags (with key ring)

Catalog Number	Description	Size
775-130	Silver	2 3/8" x 1 1/4"
775-131	Gold	2 3/8" x 1 1/4"
775-132	Red	2 3/8" x 1 1/4"
775-133	Blue	2 3/8" x 1 1/4"
775-134	Green	2 3/8" x 1 1/4"
775-135	Black	2 3/8" x 1 1/4"

Each tag allows for 5 lines of custom text max, 32 characters each.

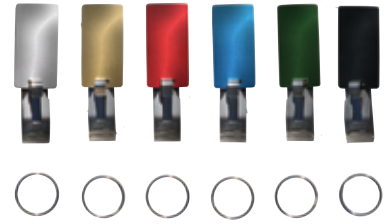


Tamper Evident Locks

Catalog Number	Description
775-400	Single use, blue, box of 1000



Indicator Cards

Catalog Number	Description
775-300	Indicates for use in gas or steam, box of 250

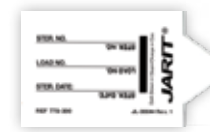


Filters

Catalog Number	Description
775-200	Single use, round 7 1/2", box of 1000



Retention Plate

Catalog Number	Description
775-212	Retention Plate

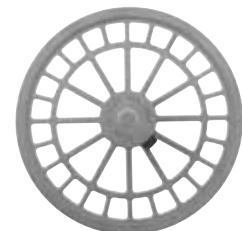


Instrument Rack for the Sterilization Container System

Features and Benefits

- 9 slot
- 5-10mm diameter silicone bars
- Compatible with the 8" Full - Size Container.
- The Rack fits inside 2" or 3" wire or perforated basket to allow ample storage underneath for other instruments and auxiliary items.
- Latex-Free and constructed with medical grade silicone grips and stainless steel frame.

Catalog Number	Description
790-010	Instrument Rack

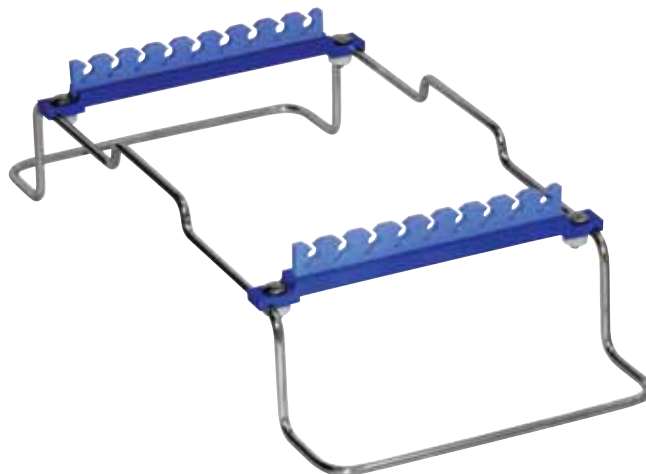

Tolerances in product specifications, due to the hand craftsmanship of the instruments, are subject to minor variations in dimensions and finishes.

For more information or to place an order, please contact:

Integra ■ 589 Davies Drive, York, PA 17402
877-468-5572 USA ■ 914-789-7094 outside USA
609-750-4259 fax ■ 914-592-8056 outside USA fax
integralife.com

